

This Page Is Inserted by IFW Operations
and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representation of
The original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

THIS PAGE BLANK (USPTO)

12. Report on the Practical Use of a New Device (Y2)**Showa 60-31706**

51. Int. Cl.4 Identification Symbol Office Reference Number

A 61 F 2/44

6779-4C

24. 44. Public Notice September 21, 1985 (Showa 60)

(5 pages)

54. Name of Idea CERAMIC SPINAL CORD PROSTHETIC MATERIAL

21. Patent Application Showa 54-162848

65. Public Exhibition Showa 56-80121

22. Application Date November 22, 1979 (Showa 54)

43. June 29, 1981 (Showa 56)

72. Inventor KITAMURA, Yasuhiro 52-11 Inoue-cho, Higashino, Yamashina-ku, Kyoto City
(c/o Kyoto Ceramic Company)72. Inventor HIRABAYASHI, Masaya 52-11 Inoue-cho, Higashino, Yamashina-ku, Kyoto City
(c/o Kyoto Ceramic Company)

71. Applicant KYOCERA 5-22 Kita Inoue-cho, Higashino, Yamashina-ku, Kyoto City

74. Agent MATSUNO, Hidehiko (Patent Attorney)

Examiner KONDO, Kanetoshi

(p.115, column 1)

57. Area covered by the Request to Register a Practical Use of a New Device

1. Made of ceramic 5, which is formed in similar shape and measurement to the spinal column intervertebral disc 7 and/or vertebral body 100; a ceramic spinal column prosthetic material that resembles a honeycomb, composed of many pillar-like holes 4 with effective openings of diameter 0.2 mm or more to allow for increase, entry and development of this bone structure 7 and/or 100 with respect to that which is inside the top and bottom surfaces 1 and 2 that come in contact with this spinal column bone structure 7 and/or 100.

2. Material that is mentioned in paragraph 1 of the "Request to Register a Practical Use of a New Device," where the above mentioned pillar-like holes 4 are the blind holes 41 that extend and stop at the inside of each of the top and bottom surfaces 1 and 2 of the above mentioned ceramic 5.
3. Part material that is mentioned in paragraph 1 of the "Request to Register a Practical Use of a New Device," where the above mentioned pillar-like holes 4 are the penetrating holes 42 that pass through the above mentioned ceramic 5 in the lengthwise direction.
4. Part material that is mentioned in paragraph 1 of the "Request to Register a Practical Use of a New Device," where the above mentioned ceramic 5 is an alumina linked crystal, alumina ceramic body, or mullite ceramic body.

Detailed Explanation

This idea involves a ceramic spinal column prosthetic material that is used in treatment in plastic surgery, and relates to spinal column prosthetic material that is made to insert in between the spinal column structure for spinal column plastic surgery.

(p.115, column 2)

With the developments in bio-engineering, treatment involving implants of material (screw pins, blades, etc.) into bone structures have become more common in the areas of dentistry and plastic surgery, taking the place of metal artificial prosthetic material; and in these cases, the applicant has provided various inventions that use ceramic, such as alumina ceramic, that, as implant material, adapt well with bone structures and are not harmful in any way. For instance, there are structures like the ceramic bone implant part material that has many holes at the contact surface area that is indicated in (1976)-116809, where it is imbedded into the bone structure to repair the bone structure through the entry of newly created bone structure through the many air passages created on the front surface area that comes in contact with the bone structure.

However, these implant material for the bone have been for hard, solid structures, and in structures that have relatively little load; thus, to use these as they are, structurally and functionally, for complex and for prosthetic material for the spinal column that has high degree of load, has been found difficult from the point of view of strength. However, the spinal column supports body weight,

and when lifting or carrying things, there is significant load placed in the lengthwise direction of the spinal column, and at the same time, because there is constant twisting force applied with the movement of the body, there has been a problem, in terms of strength, in using this type of bone implant material that involves three-dimensional space as prosthetic material.

(p.116, column 1)

For this reason, repair procedure for the spinal column has often been to insert into the patient artificial metal prosthetic material that is durable, or to remove the ilium which is of the same soft structure as the spinal column bone structure, and transplant this into the patient. However, when artificial metal prosthetic material is used, there is a problem with it not adapting biologically, and with harmful effects from corrosion within the body. Also, transplants using ilium requires time and labor with surgery, as well as cause the patient much pain and additional burden.

This idea takes into account the existing situation with the current spinal column prosthetic methods, and aims to provide an artificial prosthetic material that is superior in maintenance, has the strength to withstand sufficient load to the spinal column, has no harmful effects on the body, and adapts well to the body.

The following is a detailed explanation of the spinal column prosthetic material that is presented here.

The spinal column prosthetic material (SCPM) presented here adapts well to the body, is composed of ceramic material that has no harmful effects to the body, can be performed by simply inserting it into the appropriate section, and the point is that it is composed of many pillar-like holes in the lengthwise direction, like a honeycomb, with effective openings of diameter 0.2 mm or more to allow for increase, entry and development of this bone structure, with respect to the part that is inside the surface that comes in contact with the ceramic spinal column bone structure.

The effective openings mentioned here refers to that which is necessary for the bone structure that is adjacent to the top and bottom of the SCPM to increase and enter; it is also called the width diameter. Therefore, as long as this width diameter is of a value that allows for increase and entry of adjacent bone structures, the diameter of the pillar holes inside is not a problem.

Also, that the pillar-shaped holes have an effective opening diameter of 0.2 mm or more is a necessary requirement for the effectiveness of the SCPM, where when an adjacent bone structure enters inside the SCPM (the vertebral body that makes up the spinal column is soft compared to the harder ceramic, so the vertebral body bone structure inserts easily into the pillar-shaped spaces of the SCPM, and this fact has already been confirmed in medical experiments), and on the prosthetic area is formed compound structures due to the ceramic and bone structures that have entered and increased. Load on the spinal column in the vertical direction is significant, and there is enough strength to withstand the...

(p.116, column 2)

...force that is put in a twisting direction, and furthermore to easily conduct prosthetics that has maintenance (fixed) strength. Also, with this ceramic, because the prosthetic spinal column in general does not hinder the function it originally has, it goes without saying that it is manufactured in the well known method of "pressure-grind baking" alumina or mullite-type ceramic, with the shape that corresponds to the prosthetic part (for instance, the cervical vertebrae, thoracic vertebrae, lumbar vertebrae, and the intervertebral discs of each). Next, several ideal implementation samples of SCPM will be given, with detailed explanations.

Diagram 1 shows that which is used, among the materials for spinal column prosthetics, for prosthetic intervertebral discs. And also following, in 7 is shown a section of an intervertebral disc that has been removed, and 7 is a healthy intervertebral disc that remains inside the body even after spinal column prosthetics.

Intervertebral disc 7 ... is each of the 100 vertebral bodies that form the spinal column ... and is a fibrous soft structure that communicates to each other. From the point of view of the entire spinal column, it functions as sort of a shock absorber.

This type of incidents where prosthetic intervertebral discs become necessary are caused when, for instance, it must be removed because of problems such as a ruptured disc, and these are times when SCPM becomes especially important.

For instance, as shown in Diagram 2, simply insert the SCPM that is formed into the appropriate shape into the reciprocal spaces of vertebral body 100, 100 that has removed the intervertebral disc 7; and when only one part of vertebral body 100, 100 is removed, the suitable shape of the SCPM should be chosen to suit the part that has been removed.

Thus, when the intervertebral disc 7 undergoes prosthetics, the adjacent bone structures 100, 100 goes into ceramic 5 through holes 41 provided at the top and bottom surfaces. The SCPM 50 forms a compound structure made of ceramic and the bone structures that increase and enter by means of the principles described above, combines strongly with the adjacent top and bottom vertebral body 100, 100, and becomes a new structure that makes up a vertebral body.

In this case, in place of the intervertebral disc 7 that is a soft structure with high contractibility, a harder ceramic 5 is being used, thus it is no longer possible to fulfill a shock absorber function mentioned above, which existed in a healthy intervertebral disc. This is the same in the transplant case of the ileum (a hard structure) also mentioned above.

(p.117, column 1)

Thus, with respect to the bad effects on the body that existed at times, and poor adaptability to the living body when artificial metal prosthetic materials were used, the SCPM made from ceramic eliminates these problems; and with respect to transplanting the ileum as a means of prosthetics, there is no more need to go through the extraction surgery of the ileum that is troublesome, and a huge burden on the part of the patient. Therefore, it proves to be significantly beneficial in terms of surgery as well as physically.

Furthermore, SCPM can obviously be used for prosthetics that include those cases when the vertebral body itself is missing from the spinal column.

For instance, for prosthetics of missing areas in the spinal column, including the vertebral body, follow Diagram 3 and produce a ceramic in the shape of the desired part, and set screws 61 and 62 onto the top and bottom edges of surface 12 to which the adjacent bone structures 100, 100 are to be connected. If screws 61 and 62 are made so that they are threaded in opposite directions to one another, then screwing and securing onto the adjoining bone structure 100, 100 will be simple and sound, thus favorable.

Diagram 4 shows a prosthetic condition where the SCPM shown in Diagram 3 is put between vertebral body 100, 100.

Other ideal implementation samples of SCPM are shown in Diagrams 5, 6, 7, 8, and 9. Diagram 5 shows top and bottom surfaces 1 and 2 that combine with the bone structure, and instead of the blind hole in Diagram 1, there is ... a slot 42. Diagram 6 shows the top and bottom surfaces of SCPM shown in Diagram 1, ...with multiple rounded protrusions 3, making it easier to insert into the adjacent bone structure, and to have an intervertebral disc that has a stronger connection to the adjacent bone structure. Diagram 7 shows a prosthetic vertebral body with, ... Instead of the screws 61 and 62 at the top and bottom edges in Diagram 3, having multiple rounded protrusions 63. Diagram 8 has penetrating holes 41 whose cross-sections are shaped like that of an orange. Diagram 9 shows that with the many passage holes 41. Needless to say, these additional methods...

(p.117, column 2)

...can be applied as deemed necessary. Furthermore, if SCPM is formed to match the shape of the vertebral body part to undergo prosthetics, it is possible to perform prosthetics on various parts of the vertebral body, and especially when doing prosthetics on areas composed of the spinal column including joint areas, set an artificial joint that can freely move forward, preferably within 10°, to the middle range of the ceramic, as necessary.

As mentioned above, the SCPM is adaptable to the body structure more than the common metal prosthetics material, and the structure allows for application of the ceramic as spinal column prosthetic material that is not harmful to the body. It is effective in replacing the bone transplant since it involves simply inserting it into where spinal column prosthetics becomes necessary, and at the same time, it allows for entry, increase and development of the bone structure that is adjacent to the pillar-like holes that lies lengthwise. The ceramic and newly created bone structure form a compound structure that strongly connects to adjacent bone structures, so it maintains strong connection even under added load and lateral deviation force to the spinal column. It is an epoch-making progress in plastic surgery.

A Simple Explanation of the Diagrams

Diagram 1(a) shows one implemented sample (for intervertebral disc) using SCPM, from a diagonal view. Diagram 1(b) shows a diagonal view with one portion cut out. Diagram 2 shows the condition of the SCPM in Diagram 1(a) when used in spinal column prosthetics. Diagram 3 shows another SCPM implementation samples (for the connection part of the intervertebral disc and the vertebral body), from a diagonal view. Diagram 4 shows the condition of the SCPM in Diagram 3 when used in spinal column prosthetics. Diagram 5(a) shows another SCPM implementation sample (for intervertebral disc), from a diagonal view. Diagram 5(b) shows a diagonal view with one portion cut out. Diagrams 6, 8, and 9 are for intervertebral discs, and Diagram 7 shows the implementation sample, diagonal view, of the SCPM, showing each of the connecting parts of the intervertebral discs and the vertebral body.

Explanation of symbols 1 ... top surface, 2 ... bottom surface, 3 ... protrusion, 4 ... pillar-shaped holes, 41 ... penetrating holes, 42 ... blind holes, 5 ... ceramic, 7 ... intervertebral disc, 100 ... vertebral body, 61, 62 ... thread for screw, 63 ... protrusion.

THIS PAGE BLANK (USPTO)

⑨ 日本国特許庁(JP)

⑩ 特許出願公開

⑪ 公開特許公報(A)

昭60-31706

⑫ Int. Cl.⁴

A 47 C 7/38

識別記号

庁内整理番号

7309-3B

⑬ 公開 昭和60年(1985)2月18日

審査請求 有 発明の数 1 (全3頁)

⑭ 発明の名称 ヘッドレスト装置

⑮ 特 願 昭58-139822

⑯ 出 願 昭58(1983)7月30日

⑰ 発 明 者 齊 藤 正 次 大和市桜森3の5の10 第2コーポ福本101
⑱ 出 願 人 池田物産株式会社 横浜市保土ヶ谷区東川島町1番地の3
⑲ 代 理 人 弁理士 秋 山 修

明 細 書

1. 発明の名称

ヘッドレスト装置

2. 特許請求の範囲

シートバックの上部に、両側部のヘッドレストアームを介してヘッドレストを回動自在に設けるとともに、シートバック側の部材とヘッドレスト側の部材とに、ヘッドレストを固定するためのロック機構を設けたことを特徴とするヘッドレスト装置。

3. 発明の詳細な説明

本発明は、ヘッドレスト装置にかかり、さらに詳しくはシートバックにヘッドレストを回動自在に設けたヘッドレスト装置に関する。

従来、シートバックの上部に、回動自在にヘッドレストを取り付けられたものが提案されているが、ヘッドレストを前方に回動させ、シートバックの上部に定置した使用状態と、シートバックの後方に回動させ、シートバックの背面に嵌め込む収納状態とにセットし得るように構成されているだけ

である。したがって、前記従来の回動式のヘッドレストは、任意の回転角に調整できず、多様な用途に使用できない欠点がある。

本発明の目的は、ヘッドレストが回動可能であつて、かつ任意の回転角に調整し、固定し得るヘッドレスト装置を提供するもので、以下図面を参照して詳述する。

第1図、第2図および第3図は本発明の一実施例を示し、第4図および第5図は色々な使用形態を示す。

その第1図に示すように、シートクッション1にシートバック2がデパイヌ3を介して傾倒角を調整可能に取り付けられている。

前記シートバック2の上部には、両側部にヒンジピン4が設けられ、このヒンジピン4にヘッドレストアーム(以下、アームという)5が屈支されており、このアーム5を介してヘッドレスト6が回動自在に設けられている。尚ヘッドレスト6の使用状態における背面側にはパッド材を設けなくてもよい。

そして、シートバック2側の部材であるシートバックフレーム7と、ヘッドレスト6側の部材である一方のアーム5間に、ヘッドレスト6を任意の回転角に調整しかつ固定するためのロック機構8が設けられている。このロック機構8は、第2図および第3図に示す実施例のものは、シートバックフレーム7側に固定のラチェット9と、アーム5の根元5'側に回転可能に取り付けられかつラチェット9に掛け外し可能に設けられたフース11とを備えている。前記ラチェット9は、シートバックフレーム7に固定された非円形のラチェットビン10に嵌合され、かつ固定されている。前記フース11は、枢支ピン12を介してアーム5に回転可能に取り付けられ、かつばね13によりラチェット9にかみ合う方向に付勢されている。また、フース11にはかみ合い解除用のビン14が設けられ、このかみ合い解除用のビン14は、アーム5に枢支ピン12を中心として円弧形に形成された長穴15に挿入され、この長穴15から突出されたビン14の頭部にはノブ16が取り付けられていて、このノブ16を

操作することにより、フース11を前記ばね13に抗して第2図の矢印a方向に回転させることによりラチェット9からフース11を外し得るようになっている。

なお、前記シートバック2に固定されたラチェットビン10には、一端部から内部に向かってビン穴10'が形成され、このビン穴10'に一方のアーム5の外部からヒンジピン4が差し込まれており、ロック機構8を取り付けている側のアーム5は、このラチェットビン10に差し込まれたヒンジピン4を支点として回転し得るようになっている。また、他方のアーム5のヒンジピン4は、シートバックフレーム7に直接明けられたビン穴(図示せず)、またはシートバックフレーム7に固定されたビン(図示せず)に明けられたビン穴に当該アーム5の外部から差し込まれている。

前記実施例のヘッドレスト装置では、ヘッドレスト6を使用しない場合には第1図に示すように、シートバック2の両側部に位置するアーム5を介してシートバック2の背面側の両側部に回転させ

て収納しておく。

ついで、ヘッドレスト6を使用する場合には、ロック機構8のノブ16を操作し、フース11を第2図の矢印a方向に回転させ、ラチェット9からフース11を解除させたりして、第4図に示すように、アーム5を介してヘッドレスト6を前方に回転させ、ヘッドレスト6を使用者の体型や好みに合わせて調整し、前記ロック機構8のノブ16を放し、ばね13によりフース11をラチェット9にかみ合う方向に回転させ、ラチェット9とフース11のかみ合いを通じてヘッドレスト6を前記調整された位置に固定し、使用する。

また、デバイス3を介してシートバック2の傾倒角を変えた場合には、その傾倒角に合うように、ヘッドレスト6を調整して使用することもできる。

ついで、本発明を野外取り出しシートに適用する場合には、第5図に示すように、比較的小さな2つのシート17、18を主シート(図示せず)より取出し、この2つのシート17、18を前後に並べて用い、前方のシート17のヘッドレスト6を収納さ

せた状態でシートバック2を水平近くまで後倒させ、後方のシート18のシートバック2を任意の傾倒角に後倒させ、このシート18のヘッドレスト6を起立状態に調整して固定することによりシート18のシートバック2の実質的な高さを高くすることができ、後方のシート18のシートバック2とヘッドレスト6により使用者の上半身を支え、脚部を前方のシート17のシートバック2からシートクッション1の方に伸ばし、リラックスし得る簡易ベッドとして使用することができる。

なお、本発明においてロック機構8の具体的構造は、図面に示す実施例に限らず、要はアーム5を介してヘッドレスト6を任意の回転角に調整しかつ固定し得る構造であればよい。

以上が本発明にかかるヘッドレスト装置の実施例の構成であるが、かかる構成によればシートバックの上部に、両端部のアームを介してヘッドレストを回転自在に設けるとともに、シートバック側の部材とヘッドレスト側の部材間に、ヘッドレストを任意の回転角に調整しかつ固定するため

図中、2…シートバック、4…ヒンジピン、
5…アーム、6…ヘッドレスト、7…シートバック
フレーム、8…ロック機構。

特許出願人 他田物産株式会社
代理人 弁理士 秋 山




のロック機構を設けたことにより、ヘッドレストを任意の回転角に調整し、固定し得るので、使用者の体型や好みに合わせて調整でき、またシートバックの傾斜角を変えた場合には、その角度に適合するようにヘッドレストの回転角を調整でき、あるいは野外取り出しシートの場合には、比較的小さな2つのシートを並べて用い、一方のシートのシートバックの実質的な高さを高くする事ができ、または回転角を任意に調整し得るヘッドレスト付きの簡易ベッドとして使用可能となる等、シートとの組み合わせにおいて、色々な形態に調整し、使用し得る効果がある。

なお、上述においては本発明の一実施例を示したに留まり、本発明の精神を脱することなしに種々の変形、変更をなし得ること明らかであろう。


4. 図面の簡単な説明

第1図は本発明の一実施例を示す斜視図、第2図はロック機構の一例を示す分解斜視図、第3図は要部の拡大縦断側面図、第4図および第5図は色々な使用形態を示す斜視図である。


第1図

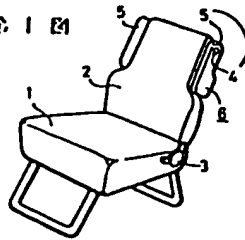


第4図

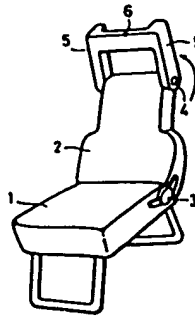


第3図

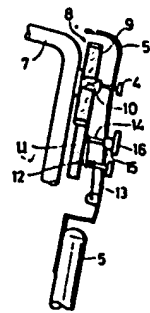

第2図

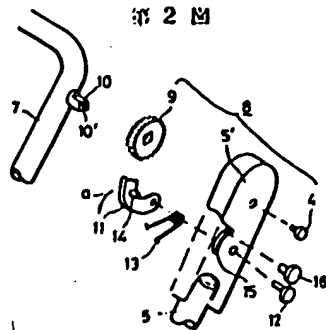

第5図

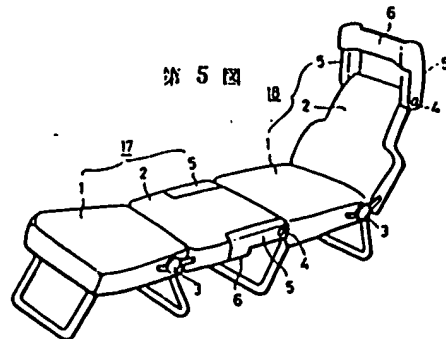

THIS PAGE BLANK (USPTO)